

TECHNICAL SPECIFICATIONS

Input Voltage	100/120/240V~, 670 VA, 50/60 Hz
Output RF Frequency	1 MHz
Fuse	T10A, 250V
Max. RF Output Power	Up to 150 watts
(MP)² Applicators	DiamondPolar™ Applicator with Integrated Thermometer OctiPolar™ Applicator with Integrated Thermometer 4D Face Applicator with Integrated Thermometer (4 suction levels) 4D Body Applicator with Integrated Thermometer (4 suction levels)
Output Voltage	Up to 130V @ Rated Load
Magnetic Pulse Frequency	15 Hz
Magnetic Field	Up to 15 Gauss
Dimensions	40 cm x 40 cm x 100 cm (D x W x H)
Weight	40 kg / 88.2 lb



الشركة الطبية العالمية
Global Medical Co.



www.gmedco.com

+ 974 4442 1661

+ 974 44429588

Building 69 , Al Juwan Street, Al Sadd,
Doha, Qatar

FACE & BODY SOLUTIONS

FOR CHALLENGING AESTHETICS



VENUS LEGACY

The highly advanced non-invasive medical device that most effectively resolves challenging face and body aesthetic needs, such as non-surgical cellulite reduction and wrinkle reduction. Treatments are comfortable with no downtime and the technology used is safe for all skin types. Venus Legacy™ is powered by proprietary (MP)² technology and features the cutting-edge VariPulse™ technology, which enhances the lipolysis effects of RF energy, increases blood circulation, and stimulates lymphatic drainage.

Zero Consumables

Zero Discomfort

Zero Downtime

TECHNOLOGY

1 | MULTI-POLAR RADIO FREQUENCY


Multi-Polar Radio Frequency (RF) uses a unique algorithm to deliver more homogeneous heating to multiple tissue depths, allowing for faster buildup of heat and easier maintenance of the therapeutic temperature. This results in greater efficacy and enhanced patient safety. With homogeneous heating, there is also no need for topical cooling agents.

2 | PULSED ELECTRO MAGNETIC FIELDS

The effect of Multi-Polar RF is enhanced by Pulsed Electro Magnetic Fields (PEMF), a non-thermal technology emitted through the applicator's electrodes. PEMF has been proven effective in triggering regenerative processes in the skin and has been used in conventional medicine for decades.


THE SYNERGISTIC EFFECT OF MULTI-POLAR RF AND PEMF

(MP)² synergistically increases collagen production through two independent mechanisms, Multi-Polar RF (thermal) and PEMF (non-thermal). The RF directly stimulates fibroblasts, while PEMF induces fibroblast proliferation through the release of the growth factor FGF-2. Together, they work to increase collagen synthesis. (MP)² is proven and effective in remodeling of collagen for skin tightening and the creation of new capillaries, which renews blood supply.

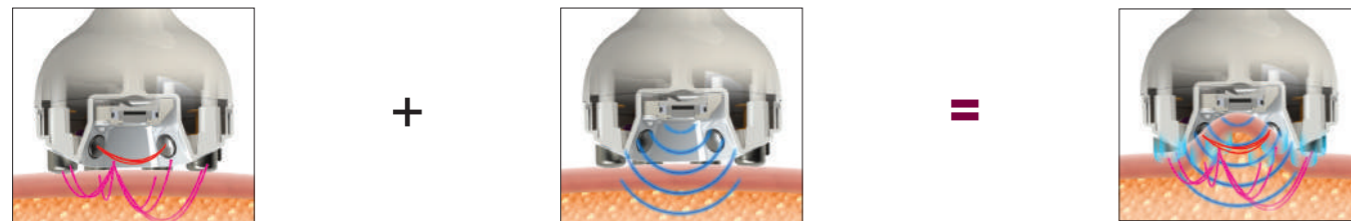

3 | VARIPULSE™ TECHNOLOGY

VariPulse™ is an adjustable pulsed suction that facilitates deep energy penetration. It enhances the penetration and effects of RF energy by improving the contact surface between electrodes and tissue, thereby improving RF coupling and heat delivery. VariPulse™ also functions to increase blood flow and enhance lymphatic drainage.

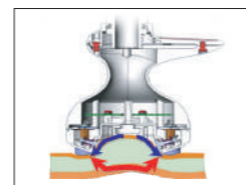

4 | REAL-TIME THERMAL FEEDBACK

All four Venus Legacy™ applicators are equipped with integrated sensors that provide real-time thermal feedback, allowing for easy, immediate, and continuous monitoring of the skin's temperature profile.



APPLICATORS

1. DIAMONDPOLAR™ APPLICATOR

The DiamondPolar applicator delivers Multi Polar Radio Frequency and Pulsed Electro Magnetic Fields to produce a dense, uniform heat matrix, stimulating the body's natural healing response. Technology increases collagen synthesis with minimized pain, resulting in progressing improvement to the skin.

2. OCTIPOLAR™ APPLICATOR

The Octipolar applicator is designed to stimulate the contraction of elastin fibers, giving skin a tighter and firmer texture. The treatment helps activate fibroblasts as well to repair damaged cells and produce healthier ones. This reduces the sagging of skin on treated areas.

3. 4D BODY APPLICATOR

The 4D Body applicator is designed with a multi-dimensional RF and PEMF energy matrix that can target specific treatment areas from multiple directions and positions to effectively stimulate fat breakdown and induce cellulite reduction. It has an innovative dual multi-polar electrode crown design, where the outer flat crown constantly controls the electro-thermal volumetric heating, while the inner deep tetra-polar crown provides intermittent vacuum-assisted electro-thermal selective heating. The unique array of eight electrodes with up to 7.5 cm distance between them provides a large spot size that allows for faster treatments. This distance, combined with the VariPulse™ technology, delivers a significant depth of penetration of up to 4.5 cm for a greater effect on fat, therefore delivering the best possible results.

4. 4D FACE APPLICATOR

The 4D Face applicator possesses all the unique technological features and benefits of the 4D Body applicator, but in a smaller configuration. It is designed to treat smaller areas of the body and enables operators to easily access hard-to-reach areas, such as around the cheekbones and neck.

